

# Nuclear extraction and fractionation of chromatin-associated proteins

The eukaryotic nucleus contains distinct biochemical environments that regulate access to the genetic information and modulate protein interactions. These include the nucleosol, euchromatin, heterochromatin, and insoluble structures such as nuclear bodies, the nucleolus, and the nuclear matrix.

Many nuclear proteins associate with DNA through charge-based interactions. Soluble nucleosolic proteins, including transcription factors and chromatin regulators, are already released under low salt conditions. More tightly bound proteins can be extracted in a single step at intermediate salt concentrations (total chromatin), or fractionated sequentially with increasing salt to distinguish between euchromatin and heterochromatin. Compacted chromatin and nucleosome-associated proteins that remain insoluble after high salt extraction become accessible through partial or complete micrococcal nuclease (MNase) digestion. The remaining pellet can be further processed by sonication or acidic extraction [SOP0026](#).

The resulting fractions can be used for immunoblotting, co-immunoprecipitation, nucleosome-dipping, or proteomics. This protocol has limited use for profiling transient or dynamic protein–chromatin interactions, which may be disrupted by salt or nuclease treatment.

Begin with isolated nuclei [SOP0001](#).

## Risk assessment

▷ Wear gloves, safety glasses, lab coat



Reviewed: Dec 10, 2023

## Procedures

### Choosing a nuclear extraction and fractionation workflow

- (1.) Choose a strategy for extracting target proteins from subnuclear compartments. If unsure, start with salt fractionation after partial MNase digestion or total chromatin extraction. [🔗](#)

| Target compartment | Accessibility | Buffer composition | Extraction strategy                     | Example markers |
|--------------------|---------------|--------------------|-----------------------------------------|-----------------|
| Nucleosol          | High          | <150 mM NaCl       | Low salt extraction                     | Sp1, Pol II     |
| Euchromatin        | Intermediate  | 150–300 mM NaCl    | Partial MNase digestion                 | H3K4me3, BRD4   |
| Heterochromatin    | Low           | >400 mM NaCl       | Partial MNase digestion                 | HP1, H3K9me3    |
| Total chromatin    | Varying       | >400 mM NaCl       | Full MNase digestion                    |                 |
| Insoluble histones | Very low      | Not applicable     | Acid extraction <a href="#">SOP0026</a> | Histone H1      |
| Nucleolar proteins | Very low      | Not applicable     | Sonication                              | NPM1            |

*Hint:* Nuclear extraction buffers use Tris-Cl by default. When planning for mass spectrometry, consider replacing with 10 mM HEPES pH 7.5 or dialyzing the extract. To enhance protein complex stability for immunoprecipitation, add 10–20% glycerol.

### >> Preparing nuclei for extraction and optional MNase digestion

- ☐ Tris-NaCl-KCl, 250  $\mu$ L [R](#)
- ☐ 0.5 M Ethylene glycol tetraacetic acid, 10  $\mu$ L [R](#)
- ☐ 200 U/mL Micrococcal nuclease, 50  $\mu$ L [R](#)

- (1.) Start with a nuclei pellet isolated from  $1.0 \times 10^6$ – $1.0 \times 10^7$  cells. Bring to  $1.0 \times 10^7$  /mL in TNK buffer. Gently flick the tube to disperse remaining clumps.

*Hint:* For continuous cell lines such as HeLa, U2OS, or 293T,  $1.0 \times 10^6$  cells typically yield sufficient material. Primary or finite cell strains may require  $1.0 \times 10^7$  cells per condition.

*Hint:* Supplement 10% glycerol to stabilize protein complexes for IP. Sucrose offers but weak protein stabilization.

*Critical:* The nuclear membrane is fragile! Gentle pipetting is sufficient to disrupt most of the pellet. Clumping after resuspension indicates lysis — discard and repeat. If iso-osmotic sucrose buffer containing Tergitol-type NP-40 (CAS 9016-45-9) was used to prepare nuclei, highly mobile nuclear proteins like HMGB1 may already be lost. [🔗](#)

- (2.) *Optional:* Supplement the TNK buffer for chromatin digestion with calcium and MNase. Incubate with gentle mixing at 37 °C, and inactivate on ice. See table for common conditions. 

| Target fraction | MNase    | Additives              | Treatment       | Inactivation |
|-----------------|----------|------------------------|-----------------|--------------|
| Total chromatin | 5.0 U/mL | 1 mM CaCl <sub>2</sub> | 20 min at 37 °C | 5 mM EGTA    |
| Euchromatin     | 1.2 U/mL | 1 mM CaCl <sub>2</sub> | 10 min at 37 °C | 5 mM EGTA    |

*Quality assurance:* Optimize MNase concentration and time for each cell type. Ideal digestion yields mostly mononucleosomes with minor di-/tri-nucleosome signal. Excessive MNase treatment will chop all nucleosomal DNA down to 10 bp fragments. 

*This is why:* EGTA preferentially chelates Ca<sup>2+</sup> leaving Mg<sup>2+</sup> in solution, a mediator of many protein-protein interactions.

- (3.) Pellet the nuclei at 400 × g for 10 min at 4 °C. If MNase was used, keep the supernatant for analysis.

*Hint:* In untreated samples, the supernatant may contain loosely associated proteins from nuclei with ruptured membrane. Discard unless required for analysis.

- (4.) Wash the pellet once with 2–5 vol TNK buffer and proceed to extraction.

 [TH12]

## >> Nuclear extraction and fractionation

- ☐ Nuclear extraction buffer, low-salt (LS), 10 mL (R) ☐ Nuclear extraction buffer, high-salt (HS), 10 mL (R)

- (1.) Prepare buffers with defined NaCl concentrations by mixing LS and HS nuclear extraction buffers:

| Nuclear extraction buffer | 0 mM NaCl | 75 mM NaCl | 150 mM NaCl | 300 mM NaCl | 400 mM NaCl |
|---------------------------|-----------|------------|-------------|-------------|-------------|
| Low-salt                  | 1.5 mL    | 0.2 mL     | 0.4 mL      | 0.8 mL      | 1.0 mL      |
| High-salt                 | 0.0 mL    | 1.4 mL     | 1.2 mL      | 0.8 mL      | 0.5 mL      |

- (2.) Resuspend the nuclear pellet in 100 µL extraction buffer per 1.0×10<sup>7</sup> nuclei. Add fresh inhibitors.

*Critical:* Use chilled buffers. Add protease inhibitors just before use to prevent degradation. 

- (3.) Incubate on ice for 30 min. For volumes over 100 µL, rotate at 4 °C.  30 min
- (4.) Separate soluble fraction by centrifugation at 400 × g for 10 min at 4 °C. Collect the supernatant.
- (5.) Repeat the extraction using buffers with increasing salt concentration on the remaining pellet. Analyze and pool fractions as needed.
- (6.) The remaining chromatin pellet can be collected at 16 000 × g for 40 min at 4 °C and resuspended for further processing such as acidic extraction or sonication.

*This is why:* Chromatinized histones remain insoluble even at 400 mM salt unless nuclei were fully digested with MNase.

 [HAW17; KNC20]

## Analyses

- Prepare samples from each fraction for SDS-PAGE and immunoblotting. Load material corresponding to  $2.0 \times 10^5$  cells or 10–20  $\mu\text{g}$  total protein per lane.
- Probe for compartment-specific markers to assess fractionation quality: Tubulin (55 kDa; cytoplasm), transcription factor Sp1 (80 kDa; nucleosol), and/or histone H3 (15 kDa; chromatin).

*Resource:* The Nuclear Protein Database (<https://webapps.igc.ed.ac.uk/npd/>) by Bickmore and Sutherland (2002) provides an overview of the diversity of many subnuclear compartments.

- *Optional:* Subject select fractions to immunoprecipitation, mass spectrometry, or nucleosome profiling.

*Critical:* If samples contain high salt or detergent concentrations, exchange the buffer by dialysis before LC-MS/MS.

## Recipes

### Tris-NaCl-KCl (TNK), pH 7.6

| Amount            | Ingredient               |         | Stock | Final |
|-------------------|--------------------------|---------|-------|-------|
| 0.5 mL            | Tris-Cl, pH 7.4          | ◇ R0055 | 1 M   | 10 mM |
| 1.0 mL            | Potassium chloride (KCl) | ◇ R0038 | 3 M   | 60 mM |
| 150 $\mu\text{L}$ | Sodium chloride (NaCl)   | ◇ R0046 | 5 M   | 15 mM |
| To 50 mL          | Water, reagent-grade     |         |       |       |

#### Tris-NaCl-KCl (TNK)

10 mM Tris-Cl, 60 mM KCl, 15 mM NaCl, pH 7.6

Date: Sign: R0143

### Micrococcal nuclease (MNase), 200 U/mL

| Amount           | Ingredient                                                          |  | Stock        | Final                |
|------------------|---------------------------------------------------------------------|--|--------------|----------------------|
| 50 U             | Micrococcal nuclease, from <i>Staphylococcus aureus</i> [9013-53-0] |  | 100–300 U/mg | 0.2 U/ $\mu\text{L}$ |
| 20 $\mu\text{L}$ | Sodium phosphate, pH 7.0                                            |  | 1 M          | 5 mM                 |
| To 4 mL          | Water, reagent-grade                                                |  |              |                      |

Dispense into 100  $\mu\text{L}$  aliquots. Stable for 1 year at  $-20^\circ\text{C}$ . **Critical:** Nuclease activity requires 1–5 mM  $\text{Ca}^{2+}$  and less than 100 mM monovalent salt! **Note:** Multiple freeze-thaw cycles are tolerated.

#### 200 U/mL Micrococcal nuclease (MNase)

0.2 U/ $\mu\text{L}$  Micrococcal nuclease, 5 mM Sodium phosphate



Date: Sign: R0144

### Ethylene glycol tetraacetic acid (EGTA), 0.5 M

| Amount   | Ingredient             |  | Stock        | Final |
|----------|------------------------|--|--------------|-------|
| 1.90 g   | Egtazic acid [67-42-5] |  | 380.35 g/mol | 0.5 M |
| To 10 mL | Water, reagent-grade   |  |              |       |

EGTA will not dissolve well until pH 8; titrate in with NaOH.

#### 0.5 M Ethylene glycol tetraacetic acid (EGTA)



#### **WARNING**

Skin irritation; Eye irritation

Date: Sign: R0145

### Nuclear extraction buffer, low-salt (LS), pH 7.4

| Amount   | Ingredient           |           | Stock | Final |
|----------|----------------------|-----------|-------|-------|
| 100 µL   | Tris-Cl              | ◇ R0055   | 1 M   | 10 mM |
| 20 µL    | MgCl <sub>2</sub>    | ◇ R0031   | 1 M   | 2 mM  |
| 40 µL    | EGTA                 | ◇ R0145   | 0.5 M | 2 mM  |
| 150 µL   | Triton™ X-100        | □ ◇ R0057 | 10%   | 0.1%  |
| 100 µL   | Inhibitor cocktail   | □ ◇ R0090 | 100 × | 1 ×   |
| To 10 mL | Water, reagent-grade |           |       |       |

Filter sterilize. Store at room temperature. **Critical:** Replace Tris-Cl with HEPES if planning for downstream mass spectrometry! **Hint:** The addition of detergent is optional, but facilitates resuspension of the chromatin fraction. To enhance protein complex stability for downstream IP, supplement the buffer with 10–20% glycerol.

#### Nuclear extraction buffer, low-salt (LS)

10 mM Tris-Cl, 2 mM MgCl<sub>2</sub>, 2 mM EGTA,  
□ 0.1% Triton™ X-100, □ 1 × Inhibitor cocktail,  
pH 7.4

Date: Sign: R0146

### Nuclear extraction buffer, high-salt (HS), pH 7.4

| Amount   | Ingredient             |           | Stock | Final |
|----------|------------------------|-----------|-------|-------|
| 100 µL   | Tris-Cl                | ◇ R0055   | 1 M   | 10 mM |
| 2.4 mL   | Sodium chloride (NaCl) | ◇ R0046   | 5 M   | 0.6 M |
| 20 µL    | MgCl <sub>2</sub>      | ◇ R0031   | 1 M   | 2 mM  |
| 40 µL    | EGTA                   | ◇ R0145   | 0.5 M | 2 mM  |
| 150 µL   | Triton™ X-100          | □ ◇ R0057 | 10%   | 0.1%  |
| 100 µL   | Inhibitor cocktail     | □ ◇ R0090 | 100 × | 1 ×   |
| To 10 mL | Water, reagent-grade   |           |       |       |

Filter sterilize. Store at room temperature. **Critical:** Replace Tris-Cl with HEPES if planning for downstream mass spectrometry! **Hint:** The addition of detergent is optional, but facilitates resuspension of the chromatin fraction. To enhance protein complex stability for downstream IP, supplement the buffer with 10–20% glycerol.


#### Nuclear extraction buffer, high-salt (HS)

10 mM Tris-Cl, 0.6 M NaCl, 2 mM MgCl<sub>2</sub>,  
2 mM EGTA, □ 0.1% Triton™ X-100,  
□ 1 × Inhibitor cocktail, pH 7.4

Date: Sign: R0147

## Resources

### Choosing a nuclear extraction and fractionation workflow

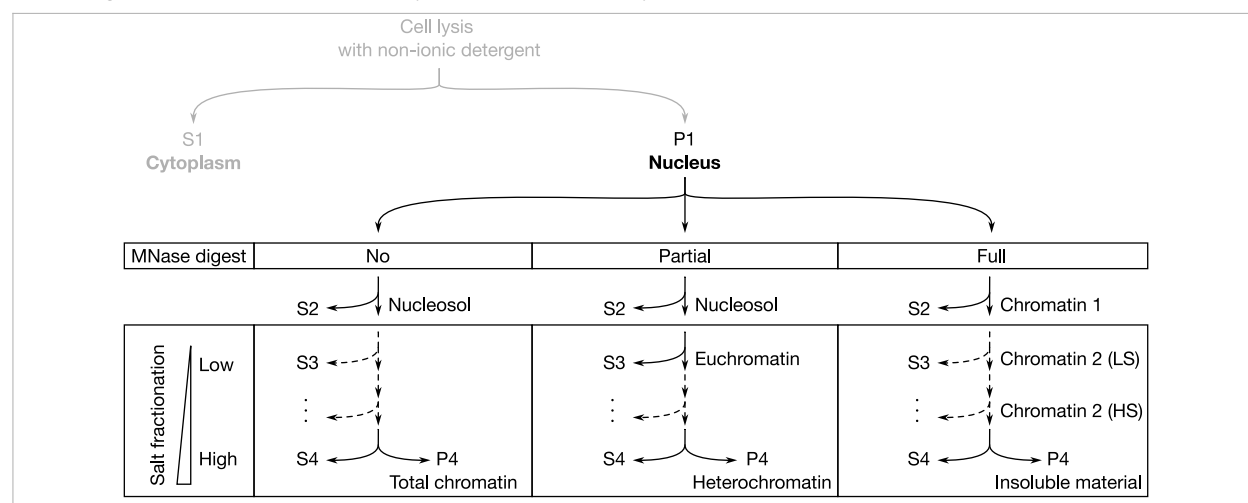

*In Step 1:* Salt extraction alone releases soluble nuclear proteins (S2). Partial MNase digestion improves solubilization of chromatin-bound proteins during salt fractionation. After complete digestion, 20% of total histone H3 and H4 may already appear in the first extract (S2), and up to 70% after the second extract (S3), regardless of salt concentration. Remaining histones in the pellet (P4) require acid extraction.

## List of references

- W.A. Bickmore and H.G.E. Sutherland, *EMBO J.* **21**(6), 1248—1254 (2002).  
 K.A. Khan, M.K. Ng, and P. Cheung, *Front. Cell Dev. Biol.* **8** 331 (2020).  
 C. Herrmann, D.C. Avgousti, and M.D. Weitzman, *Bio-protocol* **7**(6), e2175 (2017).  
 S.S. Teves and S. Henikoff, *Methods Mol. Biol.* **833** 421—432 (2012).

## Change log

- 2023-10-18 Benjamin C. Buchmuller Adaptation as SOP.  
 2023-12-10 Benjamin C. Buchmuller Revised figures; replaced washing step in isosmotic sucrose with commonly used TEK buffer; included useful hints from PMID:32457909.

From the *Lab Protocols* collection by B. C. Buchmuller and contributors (2026). Licensed under Creative Commons Attribution Share Alike 4.0 International; see <https://creativecommons.org/licenses/by-sa/4.0/> for the terms.

For research use. Provided in good faith, without warranty or liability for any use or results. Users are responsible for compliance with local regulations and institutional policies.

Current when printed. Visit <https://benjbuch.github.io/check/> or scan the QR code to check for updates.



ece3127

