

Counting cells

A hemocytometer is a specialized microscope slide with a precision-etched grid to determine the number of cells in a known volume. The most common variant for mammalian cells is the Neubauer or “Improved Neubauer” chamber. By staining with a viability dye such as Trypan blue, the proportion of live and dead cells can be determined simultaneously.

Risk assessment

- Work with human-derived material or transgenic cell lines (BSL-2)
- ▷ Wear gloves, safety glasses, lab coat
- Dispose waste after inactivation as REGULATED MEDICAL WASTE



Reviewed: Mar 14, 2023

Procedures

» Cell counting with a Neubauer cell chamber

- Neubauer cell chamber
- 0.4% Trypan blue solution (10 μ L / 1 rxn, optional)
- Hand tally

- (1.) *Optional:* Gently mix the cell suspension with an equal volume of 0.4% Trypan blue.

This is why: Trypan blue is a diazo dye excluded from live cells, allowing simultaneous viability assessment.

Hint: Automated cell counters use the same Trypan blue exclusion principle and accept 10 μ L samples in disposable slides. They reduce user-to-user variability and are preferred for high-throughput work. Follow the manufacturer’s instructions for sample preparation and calibration.

- (2.) Clean the hemocytometer and coverslip with a lint-free tissue to remove dust and debris.
- (3.) Moisten the raised glass rails of the hemocytometer. Place the coverslip over the counting surface.

Note: The coverslip is thicker than standard types to withstand surface tension and ensure correct chamber depth.

- (4.) Carefully load about 10 μ L of the sample into the chamber. The void will fill by capillary action.
- (5.) Place the hemocytometer on the microscope stage and focus. Cells should be evenly distributed and not overlapping.

- Using the hand tally, count the cells in the four corner grids of the chamber.
- Exclude cells on the left or bottom border, count cells on the right and top border (or vice versa).
- Tally unstained (live) and stained (dead) cells separately.

Quality assurance: Repeat with a more concentrated sample if the overall tally is 200 cells or less.

Hint: Most hemocytometers have nine large 100 μ m² grids that are arranged in a 3 $\times$ 3 pattern. The four corner grids with 16 small squares are used to count white blood cells and mammalian cells; the center grid is used for smaller particles such as red blood cells, platelets, or cell nuclei.

- (6.) Average the number of cells across the counted grids and multiply by 1 $\times$ 10⁴ mL⁻¹ to determine the cell concentration. Be sure to account for any dilution factor introduced by Trypan blue.

Hint: For example, if the average count across four corner grids is 150 cells and the sample was diluted 1:2 with Trypan blue, the concentration is 3.0 $\times$ 10⁶ /mL cells.

This is why: The Neubauer and “Improved Neubauer” chambers have a depth of 100 μ m between the counting grid and coverslip. Other chamber designs may differ.

Materials

Cell counting with a Neubauer cell chamber

Neubauer cell chamber, with 0.4 mm coverslip

Hand tally

Trypan blue solution ⚠ R0133 (10 µL / 1 rxn, *optional*), 0.4%

Recipes

Trypan blue solution, 0.4%

Amount	Ingredient	Stock	Final
0.4 g	Trypan Blue [72-57-1]	960.8 g/mol	0.4%
	Carcinogenicity (H350)		
10 mL	Phosphate-buffered saline (PBS), pH 7.4 ⚠ R0037	10 ×	
To 100 mL	Water, reagent-grade		

Dissolve the dye in 80% of the final volume; bring to a slow boil, cool down and make up with water.

0.4% Trypan blue solution



DANGER

Carcinogenicity

- ☐ Collect as HAZARDOUS WASTE
- ☐ DO NOT wash into sewer

Date: Sign: R0133

Troubleshooting


Cell counting with a Neubauer cell chamber

In Step 5:


- Cells are unevenly distributed under the coverslip
 - Ensure chamber and coverslip are clean and properly assembled.
 - Mix cell suspension gently but thoroughly before loading.

Resources

Cell counting with a Neubauer cell chamber

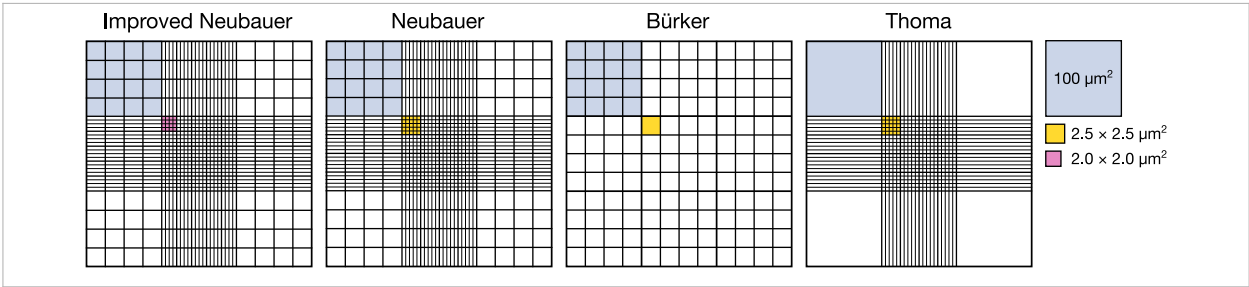

In Step 5: The design of the nine large 100 μm^2 grids varies between hemocytometer models.

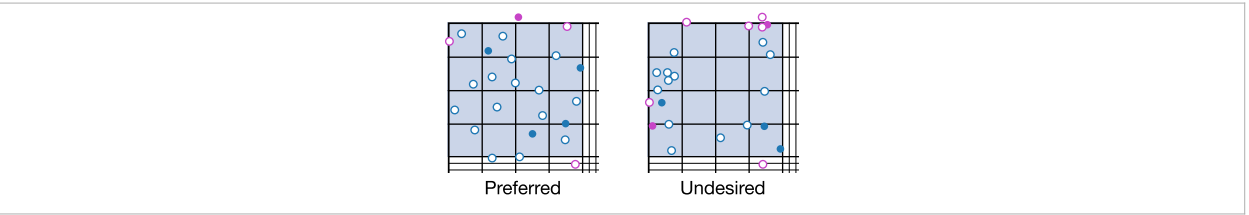

In Step 5: Particles touching the top or right borders of a grid are counted; those touching the bottom or left are excluded. Uniform cell distribution is essential.

Change log

2000-01-01	Unknown	Initial protocol; Jue Chen lab.
2015-08-07	Martin Wilson	Haemocytometer types; online resource (Agar Scientific Ltd.).
2020-05-27	Michael Haugbro	Initial protocol; Tom Muir lab.
2023-03-14	Benjamin C. Buchmuller	Adaptation as SOP.



6000308

